




EZ9Plus & EZ11Plus User Self Help:

“Resetting Atmospheric
Pressure Parameter”

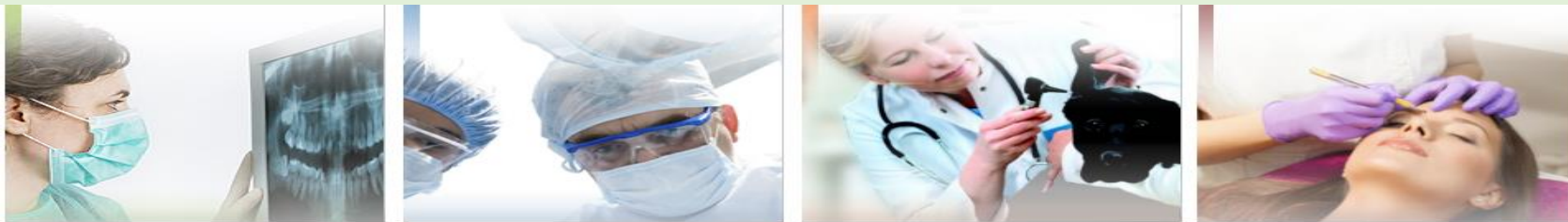


Tuttnauer[™]
Your Sterilization & Infection Control Partners



Step 1

Select any program from the Main Screen to access the Maintenance options

Step 2


Enter the programing mode by pressing the **UP** and **DOWN** keys simultaneously for 1-2 seconds and then release.





Step 1
Move the cursor to User, if it is not already highlighted, and press START/STOP key



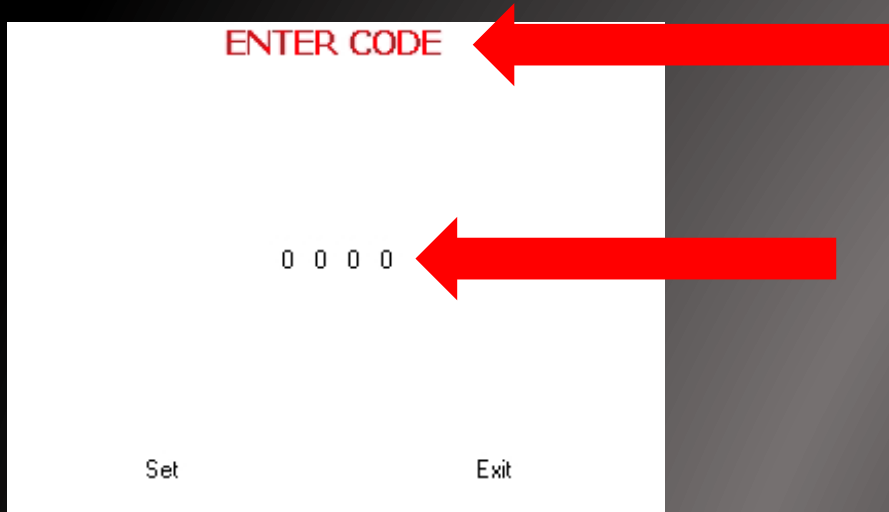

Step 2


Enter Code screen is now displayed 0000 is displayed on the screen with the cursor blinking on the right digit





Step 4


To increase or decrease the right digit, press the UP or DOWN keys

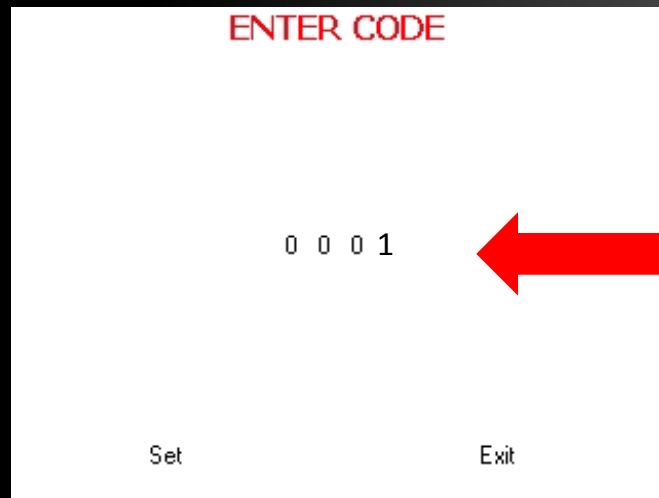

Step 5


Change the code to 0001 and move the cursor to SET by pressing the START/STOP key four times.



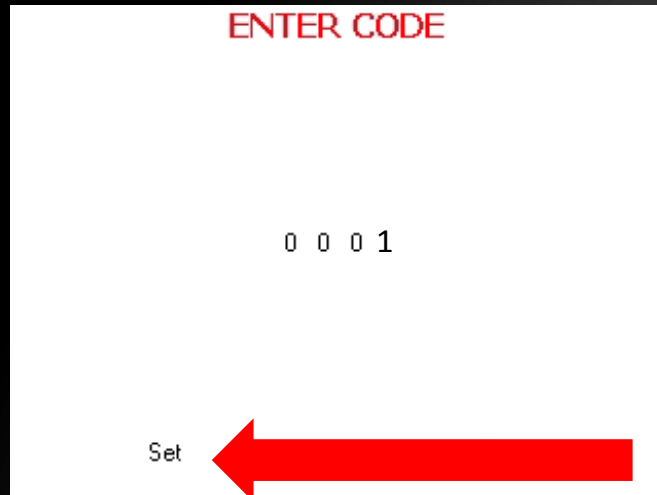

Step 6
When SET is blinking
press the UP or DOWN
key to enter the MAIN
MENU screen.



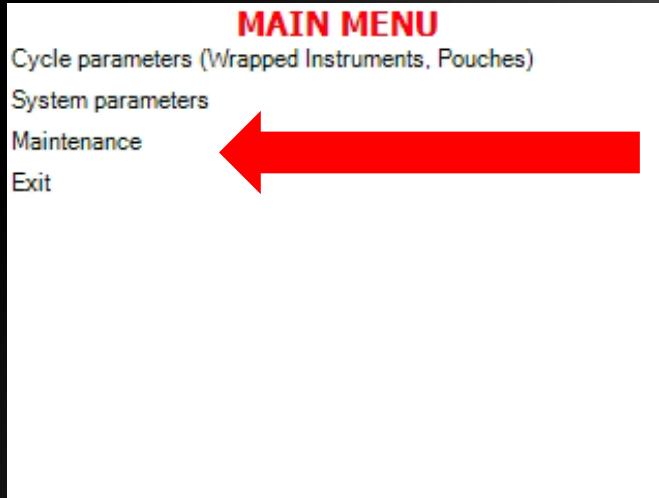


Step 7
Move the cursor to
maintenance



Step 8
press the **START/ STOP**
key.




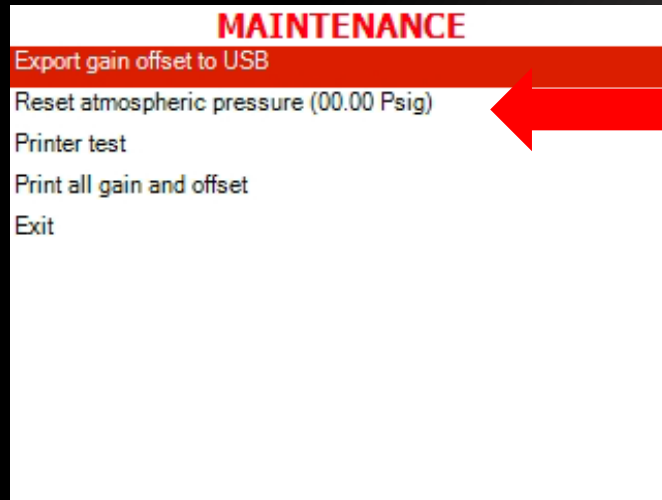


Step 9
Move the cursor to Reset atmospheric pressure then press the **START/ STOP** key.



In order to reset the atmospheric pressure make sure the door is open and the chamber temperature is less than 113°F.



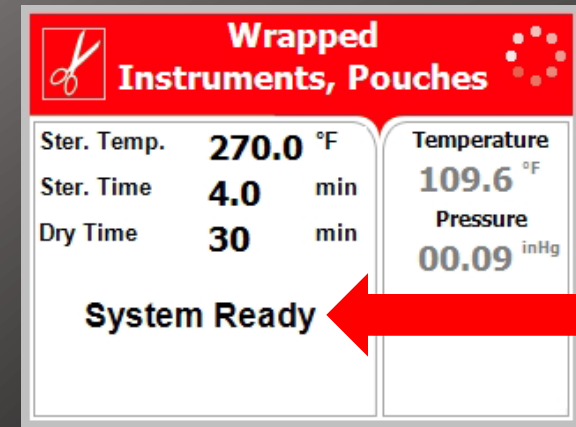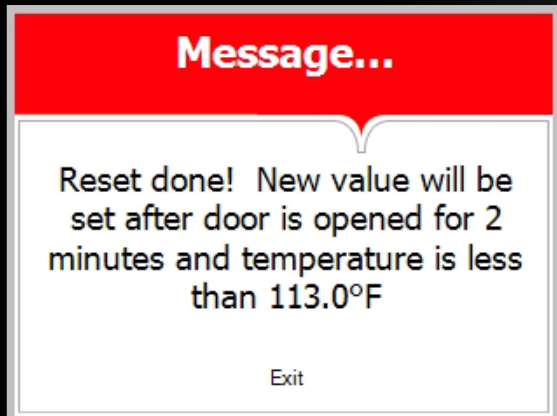

The following screen will be displayed

In order to exit this screen, press the START/STOP key

When reset is complete, the machine will return to the system ready screen

If the problem persists call your dealer or Tuttnauer USA for service.

